



Postal Address	Office
PF 060108	Kronenstrasse 19
10051 Berlin	10117 Berlin
Germany	Germany

Phone	+49-30-20 45 14 60
Fax	+49-30-20 45 14 61
E-mail	info@beromed.de
http	www.beromed.de

Office: BEROMED GMBH HOSPITAL PRODUCTS · Kronenstrasse 19 · 10117 Berlin · Germany

## Summary essential technical data - *BERO* FLOW - external catheter for male patients

**Manufacturer:** Beromed GmbH Hospital Products  
Kronenstrasse 19  
10117 Berlin  
GERMANY

**Product:** External Catheter (Type Peniflow)  
(External Catheter for urine drainage in male patients)

**Brand:** *BERO* FLOW

**UMDNS code:** 10-765

### Brief description of the product and material specification:


*BERO* FLOW is an external catheter for urine drainage in male patients suffering from temporary or chronic urinary incontinence.

Principle of operation: The catheter is placed over the tip of the penis and the sleeve is unrolled over the entire length. Adhesive strap is used around the top of the sleeve to secure the catheter to the penis. A urine drainage tube is then affixed to the catheter's distal end.

*BERO* FLOW is made of medical grade natural latex.

The product is available in sizes S (small), M (medium), L (large) and XL (X-large).

Representative drawing:

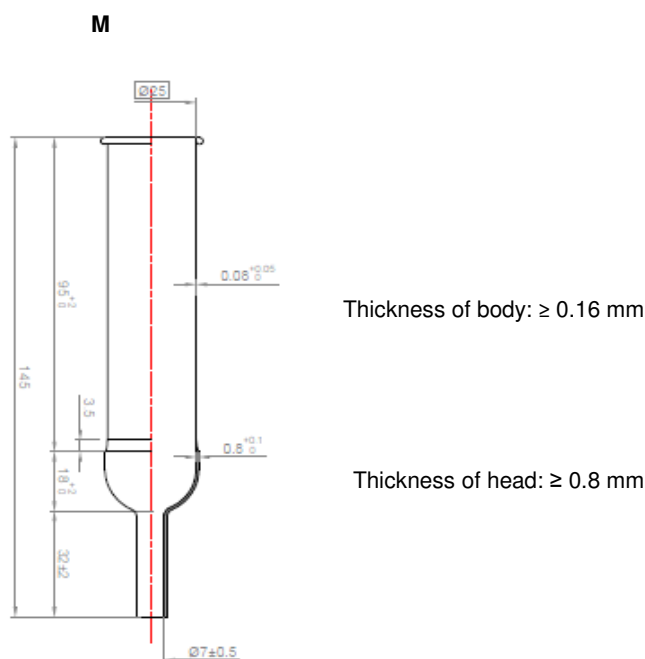

Bank	Deutsche Bank AG, Berlin
Bank Code	100 700 00
Account No.	605 37 55 00
IBAN	DE05 1007 0000 0605 3755 00
SWIFT	DEUT DE 33

Commerzbank AG, Berlin
100 800 00
905 83 45 00
DE75 1008 0000 0905 8345 00
DRES DE 33 100

Registered Office:	Berlin
Managing Director:	Jörg Bartz
Commercial Registry:	Berlin-Charlottenburg HRB 65718
VAT-No.	DE 812 404 118



<b>Postal Address</b>	<b>Office</b>
PF 060108	Kronenstrasse 19
10051 Berlin	10117 Berlin
Germany	Germany

Phone	+49-30-20 45 14 60
Fax	+49-30-20 45 14 61
E-mail	info@beromed.de
http	www.beromed.de

Office: BEROMED GMBH HOSPITAL PRODUCTS · Kronenstrasse 19 · 10117 Berlin · Germany

- Shelf life:** 5 years from the date of manufacture
- Storage:** Store in a dry, cool and dark place. Do not expose to direct sunlight.
- Disposal:** Discard the product after single use as per local medical regulations for biologically hazardous products.
- Precautions for use:** Do not use if the patient is known to develop an allergic reaction against latex.  
Do not use if the individual package is damaged.  
For single use only.  
Discard after single use to avoid risk of infection.

**Packing data:**

	pieces per
<b>Individual package</b>	1
<b>Inner box</b>	50
<b>Export carton</b>	500

**Variants:**

**BERO FLOW - External Catheter (Type Peniflow)**

Article no	Size (diameter)	
B261050	S	20 mm
B261067	M	25 mm
B261074	L	30 mm
B261081	XL	35 mm